


**Heel Incision Safety Lancets - *New and Improved Design!***


**ORIGINAL DESIGN**




**NEW AND IMPROVED DESIGN!**

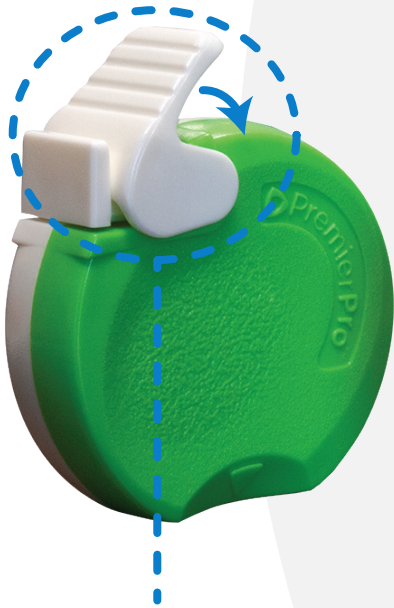

New integrated safety tab with locking trigger for one hand operation

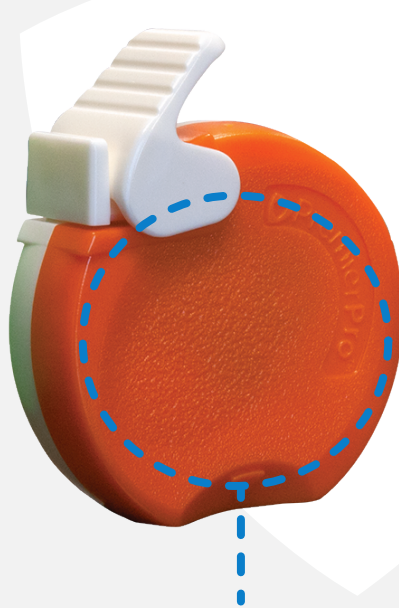

Sculpted design with anti-slip texture for easier grip

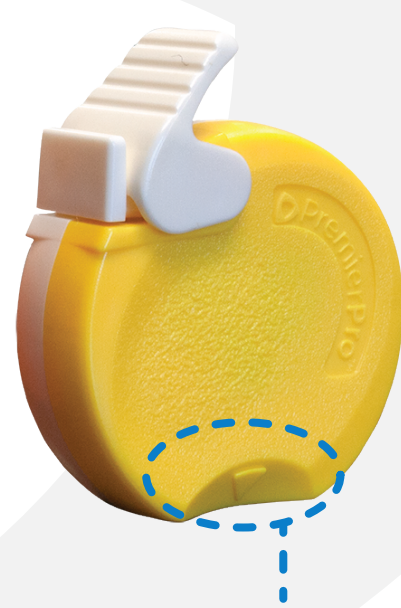

Tapered lancing position with embossed arrow to improve visibility and accuracy

**Reorder numbers and specifications have not changed:**

DESCRIPTION	COLOR	REORDER NUMBER	PACKAGING
Heel Incision Safety Lancets Depth: 0.65mm/Length: 1.4mm	Green	7511	1/pk; 50/bx; 4 bxs/cs
Heel Incision Safety Lancets Depth: 0.85mm/Length: 1.75mm	Orange	7512	1/pk; 50/bx; 4 bxs/cs
Heel Incision Safety Lancets Depth: 1mm/Length: 2.5mm	Yellow	7513	1/pk; 50/bx; 4 bxs/cs


**Pressure Activated Safety Lancets - *New and Improved Design!***
**ORIGINAL DESIGN**

**NEW AND IMPROVED DESIGN!**



Ribbed safety tab guarantees needle sterility and gives a secure grip for easy twist-off removal.

Concave housing for comfortable finger recess improves ergonomics.


No release button; gentle pressure activation.

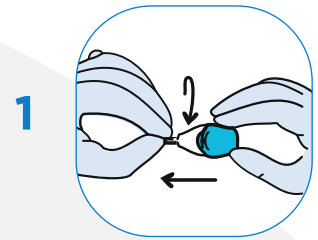

1 Twist off the protective safety tab and pull it straight out.

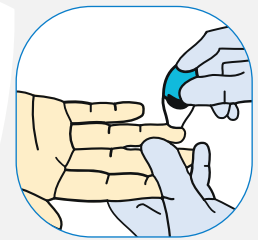

2 Press PremierPro Lancet body firmly against the puncture site to activate the device.

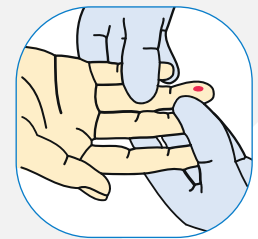

3 Gently apply intermittent pressure near the puncture site to obtain the required blood volume.


**Reorder numbers and specifications HAVE CHANGED:**
**ORIGINAL SPECIFICATIONS**

DESCRIPTION	COLOR	REORDER NUMBER	PACKAGING
Pressure Activated Safety Lancets 18G Blade x 1.8mm Depth	Green	7501	100/bx; 10 bxs/cs
Pressure Activated Safety Lancets 21G Needle x 1.8mm Depth	Light Grey	7502	100/bx; 10 bxs/cs
Pressure Activated Safety Lancets 26G Needle x 1.8mm Depth	Yellow	7503	100/bx; 10 bxs/cs
Pressure Activated Safety Lancets 30G Needle x 1.8mm Depth	Turquoise	7504	100/bx; 10 bxs/cs

**NEW SPECIFICATIONS**

DESCRIPTION	COLOR	REORDER NUMBER	PACKAGING
Pressure Activated Safety Lancets 21G Needle x 2.0mm Depth	Pink	7551	100/bx; 20 bxs/cs
Pressure Activated Safety Lancets 25G Needle x 1.8mm Depth	Blue	7552	100/bx; 20 bxs/cs
Pressure Activated Safety Lancets 30G Blade x 1.5mm Depth	Light Green	7553	100/bx; 20 bxs/cs


**Button Activated Safety Lancets - *New and Improved Design!***
**ORIGINAL DESIGN**

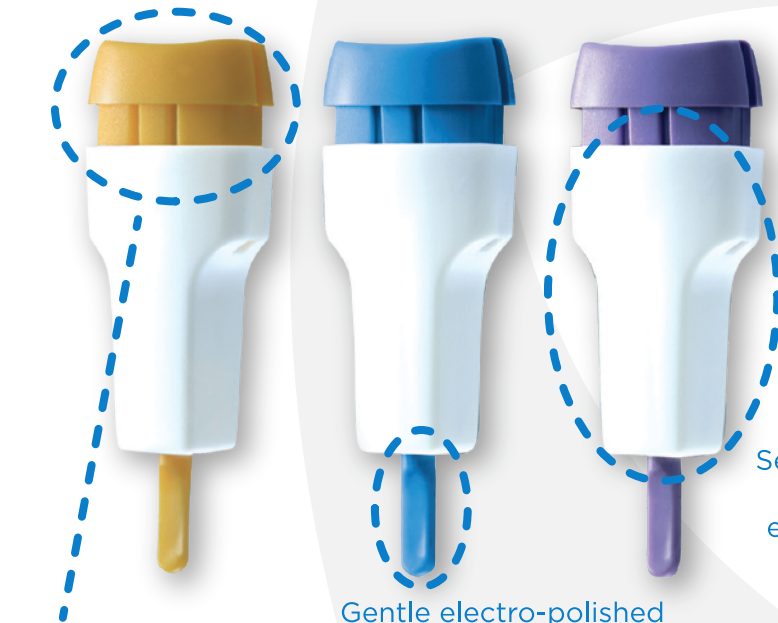**NEW AND IMPROVED DESIGN!**


Easy to handle;  
intuitive activation  
with low firing force

Gentle electro-polished  
and silicone-coated  
needle is linearly driven  
to assure lower pain  
and quick healing

Self-destructive  
mechanism  
eliminates risk  
of reuse

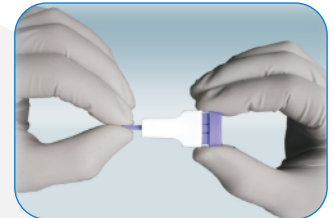**1**


Twist off the protective  
cap and pull it straight out.

**2**


Place PremierPro Lancet on  
the desired puncture site  
and press the top push  
button to activate the device.

**3**


Gently apply intermittent  
pressure near the puncture  
site to obtain the required  
blood volume.

**Reorder numbers and specifications HAVE CHANGED:**
**ORIGINAL SPECIFICATIONS**

DESCRIPTION	COLOR	REORDER NUMBER	PACKAGING
Button Activated Safety Lancets 18G Blade x 1.8mm Depth	Green	7521	100/bx; 10 bxs/cs
Button Activated Safety Lancets 21G Needle x 1.8mm Depth	Orange	7522	100/bx; 10 bxs/cs
Button Activated Safety Lancets 26G Needle x 1.8mm Depth	Yellow	7523	100/bx; 10 bxs/cs
Button Activated Safety Lancets 28G Needle x 1.8mm Depth	Purple	7524	100/bx; 10 bxs/cs

**NEW SPECIFICATIONS**

DESCRIPTION	COLOR	REORDER NUMBER	PACKAGING
Button Activated Safety Lancets 17G Needle x 2.0mm Depth	Yellow	7571	100/bx; 20 bxs/cs
Button Activated Safety Lancets 23G Needle x 1.8mm Depth	Blue	7572	100/bx; 20 bxs/cs
Button Activated Safety Lancets 28G Blade x 1.5mm Depth	Purple	7573	100/bx; 20 bxs/cs